English

For *in vitro* diagnostic use only.

Rx Only

Intended Use

The Xpert® Nasal Sample Collection Kit is designed to collect, preserve, and transport nasal swab specimens from patients with signs and symptoms of respiratory infection prior to analysis with the Xpert Xpress Flu Assay and Xpert Xpress Flu/RSV Assay.

Summary and Explanation

One of the routine procedures in the diagnosis of infections caused by respiratory viruses involves the collection and safe transportation of biological samples. This can be accomplished using the Xpert Nasal Sample Collection Kit. This system includes the Xpert Viral Transport medium that is stable at room temperature and can sustain viability (and infectivity) of respiratory viruses during transit to the testing laboratory or testing area. The formulation of the Xpert Nasal Sample Collection Kit medium includes protein for stabilization, antibiotics to minimize bacterial and fungal contamination, and a buffer to maintain a neutral pH. The Xpert Nasal Sample Collection Kit consists of one capture-cap tube containing medium and one peel pouch incorporating a sterile nylon flocked nasal specimen collection swab with a scored shaft for easy breakage. The capture-cap is designed to secure the shaft of the swab to the cap, eliminating the use of forceps to remove the swab in the laboratory or testing area.

Principles of the Procedure

The Xpert Viral Transport medium consists of modified Hank's balanced salt solution supplemented with bovine serum albumin, cysteine, gelatin, sucrose, and glutamic acid. The pH is buffered with HEPES buffer. Phenol red is used to indicate pH. Vancomycin, amphotericin B, and colistin are incorporated in the medium to inhibit growth of competing bacteria and yeast. The medium is isotonic and non-toxic to mammalian host cells. The presence of sucrose acts as a cryoprotectant, which aids in the preservation of respiratory viruses if specimens are frozen (-70°C) for prolonged storage.

Reagents

Materials Provided

Xpert Nasal Sample Collection Kit: (quantity: 100)

Each kit includes a sealed package containing:

- 1 individually wrapped sterile nylon flocked nasal swab
- 1 Xpert Viral Transport Medium tube containing:
 - Hanks' Balanced Salts
 - Bovine Serum Albumin
 - L-Cysteine
 - Gelatin
 - Sucrose
 - L-Glutamic Acid
 - HEPES Buffer
 - Vancomycin
 - Amphotericin B
 - Colistin
 - Phenol Red
 - pH $7.3 \pm 0.2 @ 25^{\circ}\text{C}$

Materials Needed but Not Provided

- Disposable gloves
- Adhesive labels or permanent markers for recording specimen information

Warnings and Precautions

- For prescription use only.
- Do not re-sterilize unused swabs.
- Do not re-pack.
- Not suitable to collect and transport microorganisms other than respiratory viruses.
- Not suitable for any other application than the intended use.
- Do not use if the swab is visibly damaged (i.e., if the swab tip is broken).
- Do not ingest the medium.
- Do not use the Xpert® Viral Transport medium for premoistening or prewetting the applicator swab prior to collecting the sample or for rinsing or irrigating the sampling sites.
- Do not bend the swab prior to specimen collection.
- Specimen collection and handling procedures require specific training and guidance.
- Observe approved biohazard precautions and aseptic techniques.
- To be used only by adequately trained and qualified personnel.
- All biological specimens should be considered potentially infectious and handled in a manner that prevents infection of laboratory personnel.
- Sterilize all biohazardous waste, including specimens, containers, and media after their use.
- Directions should be read and followed carefully

Kit Storage Requirements

The product should be stored in its original container at 2°C to 25°C . Do not overheat. Do not incubate or freeze prior to use. Improper storage will result in a loss of efficacy. Do not use after the expiration date, which is clearly printed on the outer box label

Product Deterioration


The Xpert® Nasal Sample Collection Kit should not be used if (1) there is evidence of damage or contamination to the product, (2) there is evidence of leakage, (3) the color of the medium has changed from light orange-red, (4) the expiration date has passed on any of the components, (5) the swab pouch is open, or (6) there are other signs of deterioration on any of the components.

Limitations

- Specimens should be handled aseptically.
- The Xpert Nasal Sample Collection Kit is intended for use as a collection and transport medium for respiratory viruses only.
- The Xpert Nasal Sample Collection Kit is intended to be used with the medium tube and swab provided in the kit. The use of tubes of medium or swabs from any other source could affect the performance of the product.

Specimen Collection and Handling

- Nasal specimens can be collected following the user institution's standard procedures and placed into the Xpert Viral Transport Medium (3 mL tube with transport media).
- Samples should be transported at 2–8 °C.
- Samples can be stored at 15–30 °C for up to 24 hours and at 2–8 °C for up to 7 days before processing with the Xpert Xpress Flu Assay and Xpert Xpress Flu/RSV Assay.
- Proper specimen collection, storage, and transport are critical to the performance of this test.

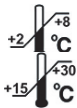Sample Collection Procedure

Proper specimen collection from the patient is extremely critical for successful isolation and identification of infectious organisms. Highest viral titers are present during the acute illness.

Note: The nylon flocked swab should not be bent prior to specimen collection.

- Collect the specimen with the swab. Take care not to touch the tip of the swab and not to lay it down.
- Aseptically remove the cap from the tube.
- Insert the swab into the tube with the transport medium.
- Break the swab shaft at the pre-scored line by bending it against the tube wall. If needed, gently rotate the Swab shaft to complete the breakage.
- Replace the cap onto the tube and close tightly.
- Label with appropriate patient information.
- Send to the laboratory or testing area for immediate analysis.

Table of Symbols

| Symbol | Meaning |
|----------------|---|
| | Catalog number |
| | CE marking – European Conformity |
| | <i>In vitro</i> diagnostic device |
| | Contains sufficient for <n> tests |
| | Sterilized using ethylene oxide |
| | Authorized representative in the European Community |
| | Do not reuse |
| | Manufacturer |
| | Temperature limitation |
| | Batch code |
| | Caution |
| Rx Only | Caution: Federal law restricts this device to sale by or on the order of a (licensed healthcare practitioner) |
| | Consult instructions for use |
| | Expiration date |
| | Biological risks |
| | Handle with care |

For Information Only - Not a Controlled Copy

Technical Support Contact

| Region | Telephone | Email |
|---|------------------------------------|------------------------------|
| US | +1 888.838.3222 | techsupport@cepheid.com |
| Australia and New Zealand | +1800 130 821 +0800 001 028 | techsupportanz@cepheid.com |
| Brazil and Latin America | + 55 11 3524 8373 | latamsupport@cepheid.com |
| China | +86 021 5406 5387 | techsupportchina@cepheid.com |
| France | +33 563 825 319 | support@cepheideurope.com |
| Germany | +49 69 710 480 480 | support@cepheideurope.com |
| India, Bangladesh, Bhutan, Nepal and Sri Lanka | +91 11 48353010 | techsupportindia@cepheid.com |
| Italy | +39 800 902 567 | support@cepheideurope.com |
| South Africa | +27 861 22 76 35 | support@cepheideurope.com |
| United Kingdom | +44 3303 332 533 | support@cepheideurope.com |
| Other European, Middle East and African countries | +33 563 825 319 +971 4 253 3218 | support@cepheideurope.com |
| Other countries not listed above | +1 408.400.8495 | techsupport@cepheid.com |

*You may also contact your local Technical Support office at: <http://www.cepheid.com/us/support/support-solutions/technical-support-offices>.

Cepheid
904 Caribbean Drive
Sunnyvale, CA 94089
USA
Phone: +1.408.541.4191
Fax: +1.408.541.4192

Cepheid Europe SAS
Vira Solelh
81470 Maurens-Scopont
France
Phone: +33 563 825 300
Fax: +33 563 825 301

For Information Only - Not A Controlled Copy

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